Hardware User's Manual

TMP812RS

Digital Thermometer

References:

TMP812RS (76-0065)

Version:

V29/06/2015

Limitation of Liability

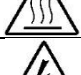
PANLAB does not accept responsibility, under any circumstances, for any harm or damage caused directly or indirectly by the incorrect interpretation of what is expressed in the pages of this manual.

Some symbols may have more than one interpretation by professionals unaccustomed to their usage.

PANLAB reserves the right to modify, in part or in total, the contents of this document without notice.

1. SYMBOLS TABLE

Recognising the symbols used in the manual will help to understand their meaning:

DESCRIPTION	SYMBOL
Warning about operations that must not be done because they can damage the equipment	
Warning about operations that must be done, otherwise the user can be exposed to a hazard.	
Protection terminal ground connection.	
Warning about a hot surface which temperature may exceed 65°C	
Warning about a metal surface that can supply electrical shock when it's touched.	
Decontamination of equipments prior to disposal at the end of their operative life	
Waste Electrical and Electronic Equipment Directive (WEEE)	

2. GOOD LABORATORY PRACTICE

Check all units periodically and after periods of storage to ensure they are still fit for purpose. Investigate all failures which may indicate a need for service or repair.

Good laboratory practice recommends that the unit be periodically serviced to ensure the unit is suitable for purpose. You must follow preventive maintenance instructions. In case equipment has to be serviced you can arrange this through your distributor. Prior to Inspection, Servicing, Repair or Return of Laboratory Equipment the unit must be cleaned and decontaminated.

Decontamination prior to equipment disposal

In use this product may have been in contact with bio hazardous materials and might therefore carry infectious material. Before disposal the unit and accessories should all be thoroughly decontaminated according to your local environmental safety laws.

3. UNPACKING AND EQUIPMENT INSTALATION

WARNING: Failure to follow the instructions in this section may cause equipment faults or injury to the user.

- A. No special equipment is required for lifting but you should consult your local regulations for safe handling and lifting of the equipment.
- B. Inspect the instrument for any signs of damage caused during transit. If any damage is discovered, do not use the instrument and report the problem to your supplier.
- C. Ensure all transport locks are removed before use. The original packing has been especially designed to protect the instrument during transportation. It is therefore recommended to keep the original carton with its foam parts and accessories box for re-use in case of future shipments. Warranty claims are void if improper packing results in damage during transport.
- D. Place the equipment on a flat surface and leave at least 10 cm of free space between the rear panel of the device and the wall. Never place the equipment in zones with vibration or direct sunlight.
- E. Once the equipment is installed in the final place, the main power switch must be easily accessible.
- F. Only use power cords that have been supplied with the equipment. In case that you have to replace them, the spare ones must have the same specs that the original ones.

- G. Make sure that the AC voltage in the electrical network is the same as the voltage selected in the equipment. **Never connect the equipment to a power outlet with voltage outside these limits.**

WARNING

For electrical safety reasons you only can connect equipment to power outlets provided with earth connections .

This equipment can be used in installations with category II over-voltage according to the General Safety Rules.

The manufacturer accepts no responsibility for improper use of the equipment or the consequences of use other than that for which it has been designed.

PC Control

Some of these instruments are designed to be controlled from a PC. To preserve the integrity of the equipment it is essential that the attached PC itself conforms to basic safety and EMC standards and is set up in accordance with the manufacturers' instructions. If in doubt consult the information that came with your PC. In common with all computer operation the following safety precautions are advised.

- WARNING**
- To reduce the chance of eye strain, set up the PC display with the correct viewing position, free from glare and with appropriate brightness and contrast settings
 - To reduce the chance of physical strain, set up the PC display, keyboard and mouse with correct ergonomic positioning, according to your local safety guidelines.

4. MAINTENANCE

WARNING: Failure to follow the instructions in this section may cause equipment fault.

- **PRESS KEYS SOFTLY** – Lightly pressing the keys is sufficient to activate them.
- Equipments do not require being disinfected, but cleaned for removing urine, faeces and odour. To do so, we recommend using a wet cloth or paper with soap (which has no strong odour). **NEVER USE ABRASIVE PRODUCTS OR DISSOLVENTS.**
- **NEVER** pour water or liquids on the equipment.
- Once you have finished using the equipment turn it off with the main switch. Clean and check the equipment so that it is in optimal condition for its next use.
- The user is only authorised to replace fuses with the specified type when necessary.

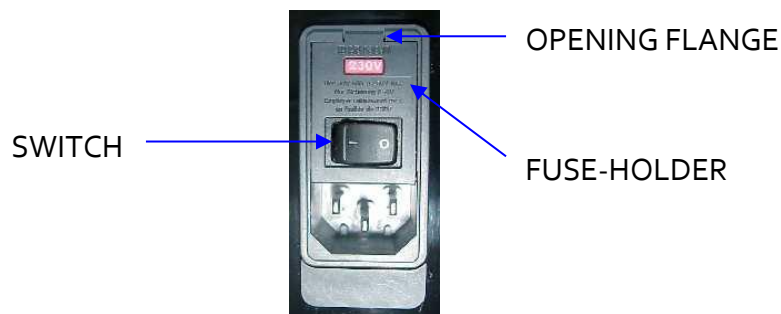
Figure 1. Power inlet, main switch and fuse holder.

FUSE REPLACEMENT OR VOLTAGE SETTING CHANGE

In case of an over-voltage or other incident in the AC net making it impossible to turn on the equipment, or if the equipment voltage setting is incorrect, check fuses according to the following procedure.

- 1 Remove power cord from the power inlet.

- Open fuse-holder by pulling the flange with a regular screwdriver.

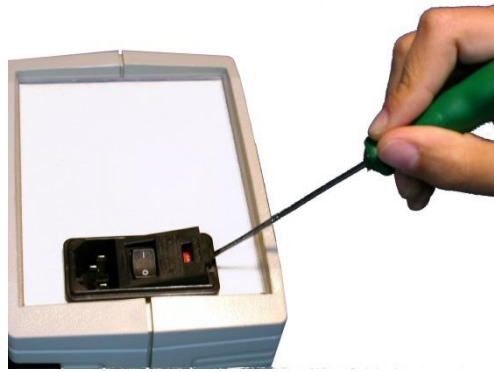
Figure 2. Open fuse-holder door.

- Extract fuse holder using the screwdriver.

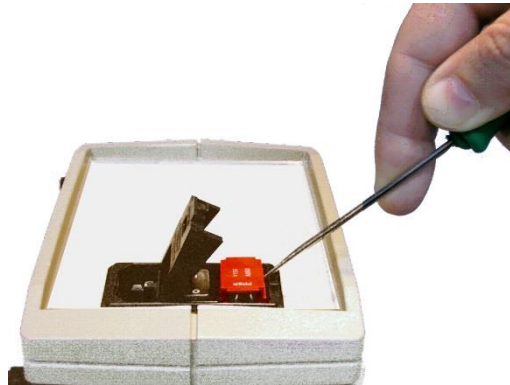
Figure 3. Extract fuse-holder.

- Replace fuses if necessary. Insert fuses in the fuse-holder in the correct position.

CORRECT

INCORRECT

Figure 4. Fuses position.

- Insert the fuse-holder again, positioning it according to the voltage in the AC net.

115V POSITON

230V POSITION

Figure 5 Fuse holder position.

- If the fuses blow again, unplug the equipment and contact technical service.

For electrical safety reasons, never open the equipment. The power supply has dangerous voltage levels.

5. TABLE OF CONTENTS

1. SYMBOLS TABLE	2
2. GOOD LABORATORY PRACTICE	2
3. UNPACKING AND EQUIPMENT INSTALATION	3
4. MAINTENANCE	5
5. TABLE OF CONTENTS	7
6. INTRODUCTION	8
7. EQUIPMENT DESCRIPTION	9
7.1. CONTROL UNIT FRONT PANEL	9
7.2. CONTROL UNIT REAR PANEL	10
8. EQUIPMENT CONNECTION	11
9. WORKING WITH THE EQUIPMENT	12
10. SENDING DATA TO A COMPUTER (SEDACOM)	14
11. TROUBLESHOOTING	15
12. PREVENTIVE MAINTENANCE	16
13. TECHNICAL SPECIFICATIONS	17

6. INTRODUCTION

The TMP 812 RS is a self-calibrated digital thermometer that can work with up to 12 temperature probes simultaneously. The RS232 output of the device can be used for the transfer of data to an external device: printer or PC associated to the SEDACOM software (both options are not included with the device, they should be purchased separately).

Figure 6. TMP 812 RS.

The analogical to digital conversion is performed with 13 bits precision.

Range of temperature readings	0°C - 50°C
Temperature precision	0.1°C
Required probe model	YSI series 400
Probe temperature reading	250 msec/probe
Renewal of the reading of each probe	3 sec.
Information format	ASCII

7. EQUIPMENT DESCRIPTION

7.1. CONTROL UNIT FRONT PANEL

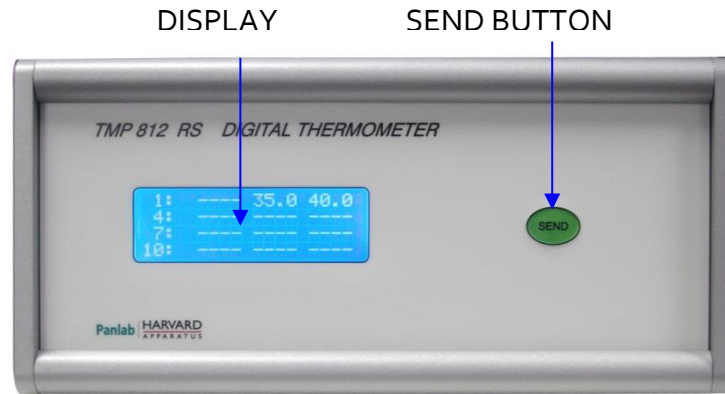
Figure 7. TMP 812 RS Front Panel.

- **DISPAY:** The display simultaneously shows the temperature of the 12 probes in °C. If a probe is broken or not connected the display will show "----" for it.
- **SEND BUTTON:** When this button is pressed, the TMP812 RS will send the temperatures to the connected external options (printer or software) through the RS2323 output.

7.2. DISPLAY

Figure 8. Digital Display.

The temperature of the 12 probes is shown on the display simultaneously in 4 rows and 3 columns. The number at the beginning of each row is the number of the first probe in that row. Thus, the first row contains probes 1, 2 and 3; the second row contains probes 4, 5 and 6; the third row contains probes 7, 8 and 9 and the last row contains probes 10, 11 and 12.

7.3. CONTROL UNIT REAR PANEL

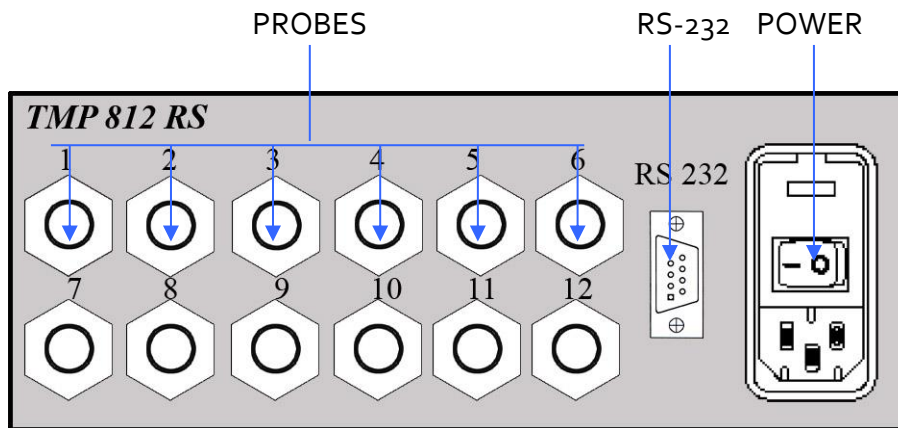
Figure 9. Rear Panel.

- **PROBES:** There are 12 jack plugs distributed in 2 rows and 6 columns. They are all labelled with probe numbers from 1 to 12. When a probe is not connected or damaged the display will show "----".
- **RS 232:** A DB9 female connector used to connect the TMP 812 RS to the computer serial port. Temperatures are sent external output device options through this connector when the **SEND** button is pressed.
- **POWER:** Power inlet, main switch and fuse holder.

8. EQUIPMENT CONNECTION

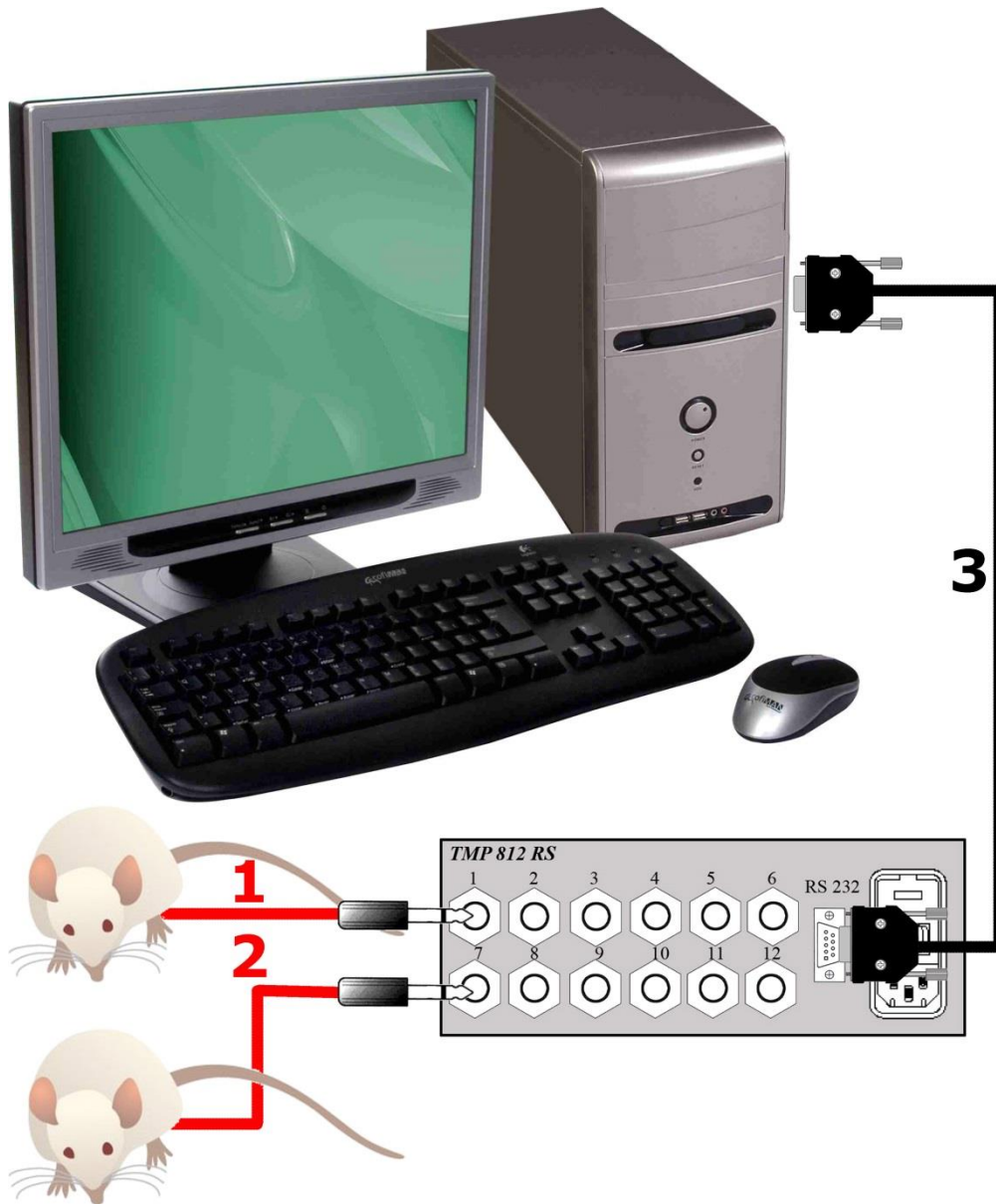
Figure 10. Connection example.

The above figure shows an example of a connection between the TMP 812 RS and the computer for temperature measurement of 2 rodents. The necessary connections and cables are listed in the following table.

	FROM	TO	CABLE
1	TMP812RS probe 1*	Rodent 1*	YSI temperature probe
2	TMP812RS probe 2*	Rodent 2*	YSI temperature probe
3	TMP812RS RS232	Computer Serial port	RS232 cable

*N is the probe number from 1 to 12.

9. WORKING WITH THE EQUIPMENT

9.1. MEASURING PROCEDURE

1. Plug in the POWER CORD, making sure to have selected the correct voltage, with the power button in its OFF position.
2. Connect the desired number of PROBES to the rear panel SOCKETS. The instrument will automatically detect which channels are working.
3. Turn on the instrument by switching the power button to the ON position. If no probe is connected, the display will show "----" for each probe. If a probe is connected, the registered temperature will appear on the display in the relevant placement.
 - a) If any probe temperature is below 0 °C, the display will show in the relevant position "<0C".
 - b) If any probe temperature is above 50 °C, the display will show in the relevant position ">50C".
4. If all the probes (12) are connected, it will take 3 seconds until every register shows its present value (it takes 250 msec. to take each register).

9.2. TESTING EQUIPMENT CALIBRATION

Along with the equipment are delivered two probes pattern of 35 ° C and 40 ° C. These probes are in fact fixed resistors that do not vary over time.

The equipment is calibrated in factory and the user cannot calibrate it. What if the user can do is to check its calibration with both pattern probes, for it simply they must be connected in any two channels and watch the display will show the values of 35 ° C and 40 ° C in the respective channels.

9.3. CLEANING THE TEMPERATURE PROBE

Read the temperature probe manual in order to follow the manufacturer recommendations.

9.4. TEMPERATURE PROBE DISINFECTION

Read the temperature probe manual in order to follow the manufacturer recommendations.

9.5. TEMPERATURE PROBE STERILIZATION

Read the temperature probe manual in order to follow the manufacturer recommendations.

10. SENDING DATA TO A COMPUTER (SEDACOM)

The purchase of the **Sedacom** software is needed for transferring the data to a computer (please contact your local provider for more information). The **Sedacom** software reference is composed by a USB Flash key containing the software Installer, a License for use and **Sedacom's** User's Manual). Follow next instructions:

- Please refer to the **Sedacom's** User's Manual for instructions on how to install and use the software with the present device.
- A serial port (RS232) communication cable (provided with the present device) is needed for connecting the present device to the computer in which the **Sedacom** software is installed. Please refer to the present User's Manual chapter 8 for instructions on how to connect this cable to the device.
- If the computer does not have any serial port, the RS232/USB adapter is needed (ref. CONRS232USB, contact your local provider for more information)

11. TROUBLESHOOTING

This table features instructions to solve the most frequent problems.

PROBLEM	SOLUTION
The equipment does not start up.	<ul style="list-style-type: none"> • Ensure that the voltage of mains is the same as that selected in the fuse holder. • Check the condition of the fuses.
The display shows "---" in the temperature of a channel.	<ul style="list-style-type: none"> • Make sure the temperature probe is properly connected on the rear panel. • If there is a probe connected, change the probe of channel to see if the problem is in the probe or in the channel.
The display shows "<0C" in the temperature of a channel.	<ul style="list-style-type: none"> • The measured temperature is lower than 0°C. • If not, change the probe of channel to see if the problem is in the probe or in the channel.
The display shows "> 50C" in the temperature of a channel.	<ul style="list-style-type: none"> • The measured temperature is higher than 50°C. • If not, change the probe of channel to see if the problem is in the probe or in the channel.
The equipment does not send data to Sedacom32.	<ul style="list-style-type: none"> • Make sure your equipment is connected to PC via RS-232. • Check that Sedacom32 settings so that the serial port and the equipment selected are correct. • Restart the equipment and the PC to do a RESET in communications.

12. PREVENTIVE MAINTENANCE

	EXPERIMENT	WEEKLY	MONTHLY	SERVICE
CHECKING CALIBRATION			<input checked="" type="checkbox"/>	
CLEANING THE TEMPERATURE PROBE	<input checked="" type="checkbox"/>			
DESINFECTING THE TEMPERATURE PROBE		<input checked="" type="checkbox"/>		
STERILIZING THE TEMPERATURE PROBE ¹				<input checked="" type="checkbox"/>

¹You must sterilize the temperature probe before sending it back for servicing.

13. TECHNICAL SPECIFICATIONS

POWER SUPPLY Input Voltage Frequency Fuses Maximum Power Conducted noise	115 /230V~ 50 /60 Hz 2 fuses 5x20mm 250mA 250V Fast 12 W EN55022 /CISPR22/CISPR16 class B
BEHAVIOURAL CONDITIONS Working temperature Working relative humidity Storage temperature	10°C a +40°C 0% a 85% RH, without condensation 0°C a +50°C, without condensation
TEMPERATURE MEASUREMENT Range of temperature readings Temperature precision Probe model required Probe temperature reading Renewal of the reading of each probe	0°C to 50°C 0.1°C YSI series 400 250 msec/probe 3 sec.
COMUNICATIONS OUTPUT Standard Interface Connector	RS232C Delta 9 contacts female connector
DIMENSIONS Width x Height x Depth Weight	232mm x 111mm x 297mm 3.5 kg

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Cornellà de Llobregat, Spain
 30/04/2014

(GB) Note on environmental protection:

After the implementation of the European Directive 2002/96/EU in the national legal system, the following applies:

Electrical and electronic devices may not be disposed of with domestic waste. Consumers are obliged by law to return electrical and electronic devices at the end of their service lives to the public collecting points set up for this purpose or point of sale. Details to this are defined by the national law of the respective country. This symbol on the product, the instruction manual or the package indicates that a product is subject to these regulations. By recycling, reusing the materials or other forms of utilising old devices, you are making an important contribution to protecting our environment.

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